

AUTOMIX

Automatic mixing machine for dental impression materials to obtain a homogeneous mixture without bubbles and precise dosing.



Why choose R&S automatic mixer ?

- + Polyvalent**
Wide choice of speed : reduced, fast and very fast.
- + Ease of use**
 - Clear labeling and automatic system (no need to keep your finger on the button).
 - Secured system: net stop of the flowing.
- + Compatibility**
 - Compatible with Turbomix cartridges and with most international brand cartridges such as 3M, Zhermack, DMG, Coltene, Walldent, Heraeus Kulzer...
 - Compatible with cartridge volume of 360 or 380 ml, ratio 5:1.
- + Ergonomic and aesthetic**
 - Small size, light and stable.
 - Smart design.




Advice from our dental experts

- >** Very stable device, allowing easy placement and removal of the impression cartridges, with clear commands, and a net stop, without flowing of the impression material.

Technical specifications and instructions for use

Dimensions (depth x width x height, in cm)	43,5 x 17 x 29	Noise level (general weighted acoustic pressure, in dB)	72
Speeds	Reduced Fast Very fast	Maximum output rating (W)	150
Weight (kg)	7	Storage	Store and transport the unit between - 20°C to 60°C (-4°F to 140°F) and a maximum relative humidity of 80%. When transporting the unit we recommend using the original packaging..
Environmental temperature (°C)	+5 to 40	Warranty	1 year
Recommended operating temperature (°C)	+18 to 40	Product reference	205-325
Power supply (V)	100/240		
Frequency (Hz)	60/50		

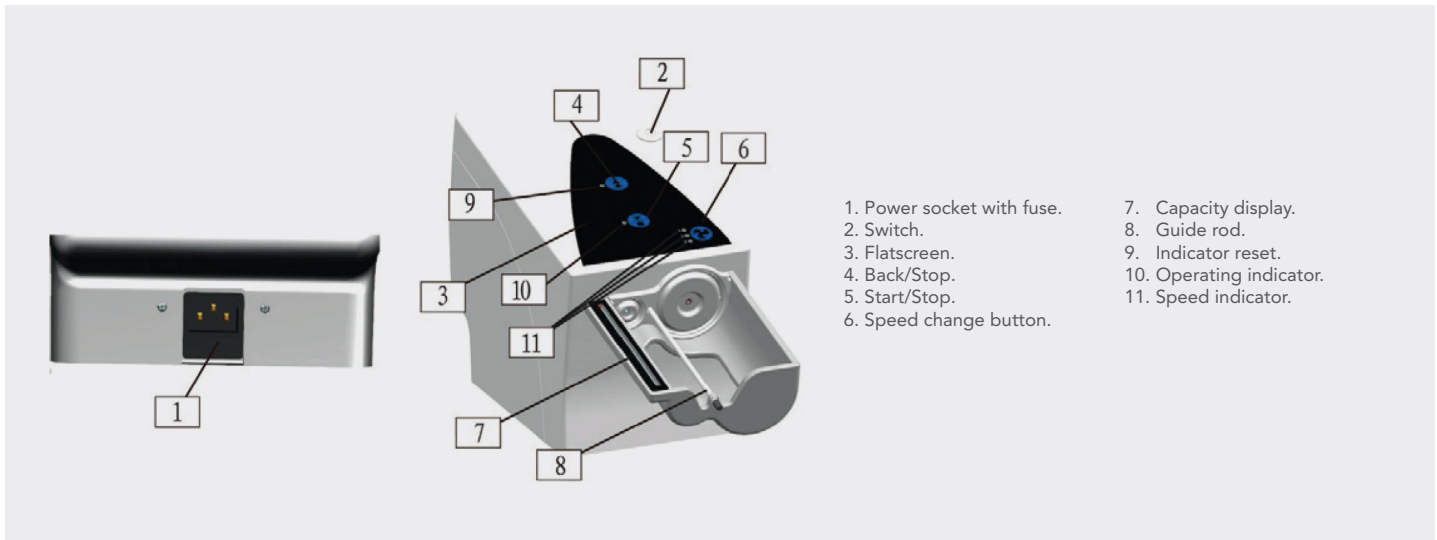

1. Power socket with fuse.
2. Switch.
3. Flatscreen.
4. Back/Stop.
5. Start/Stop.
6. Speed change button.
7. Capacity display.
8. Guide rod.
9. Indicator reset.
10. Operating indicator.
11. Speed indicator.

Complementary products



Turbomix putty (normal/fast)
Turbomix monophase fast
Turbomix heavy fast



Perforated impression trays



Turbomix mixing tips



Clamping ring